


WIRED/WIRELESS	
PROGRAMS CONFIGURATION	
CHANNELS	Number of channels of stimulation
PROGRAM CATEGORIES	TYPE OF PROGRAMS
REHABILITATION I	Atrophy, Reinforcement, Motor point...
PAIN RELIEF I	Tens Gate Control, Endorphins...
VASCULAR	Heavy legs, Capillarization...
CONDITIONING I	Resistance, Strength, Recovery...
REHABILITATION II	ACL, Rotator Cuff, Lumbar Stabilization...
PAIN RELIEF II	Low back Pain, Cervical Pain, Epicondylitis...
NEURO REHABILITATION	Spacticity, Hemiplegic Foot & Shoulder...
CONDITIONING II	Endurance, Explosive Strength, Massage...
INCONTINENCE	Stress, Urge, Mixed, Post Partum Prevention...
DENERVATED	Partial and Total, Automatic and Manual...
DIRECT CURRENT	Iontophoresis, Edema, Hyperhidrosis...
ADDITIONAL FUNCTIONS	
2+2 OR +TENS OPTION	Combine 2 different programs
MOTOR POINT PEN	Localize muscle motor point for optimal electrode placement
SKIP	Move directly to the next phase
PROGRAMMABILITY	Create your own programs
MI-SCAN	Customize the stimulation for comfort and efficiency
MI-TENS	Work efficiently behind the muscle threshold
MI-RANGE	Work efficiently at good intensity for low frequency programs
MI-ACTION	Active sessions with voluntary contraction detected by the device
MULTI-SESSIONS™	2 remotes + Smart docking station to manage several sessions independently
TRIGGER	Starts voluntary working phase
SOFTWARE	Create your patient and history lists
AUTOMATIC RAMPE SEARCH	Partial denervation
WARRANTY	3 years

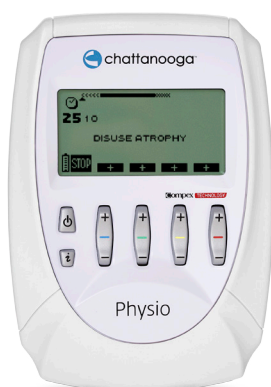

REHAB	THETA	PHYSIO	WIRELESS PRO 2CH STANDARD	WIRELESS PRO 4CH STANDARD	WIRELESS PRO 4CH FULL	
Wired	Wired	Wired	Wireless	Wireless	Wireless	
Standard	Full	Optimum	Standard	Standard	Full	
4	4	4	2 Modules	4 Modules	4 Modules	
X	X	X	X	X	X	STANDARD
X	X	X	X	X	X	
X	X	X	X	X	X	
X	X	X	X	X	X	
	X	X			X	FULL
	X	X			X	
	X	X			X	
	X	X			X	
		X				OPTIMUM
		X				
		X				
X	X	X	X	X	X	
X	X	X	X	X	X	
X	X	X	X	X	X	
X	X	X	X	X	X	
Mi-ready	2 Mi-cables	4 Mi-cables	2 Ch	4 Ch	4 Ch	
Mi-ready	2 Mi-cables	4 Mi-cables	2 Ch	4 Ch	4 Ch	
Mi-ready	2 Mi-cables	4 Mi-cables	2 Ch	4 Ch	4 Ch	
Mi-ready	2 Mi-cables	4 Mi-cables	2 Ch	4 Ch	4 Ch	
				X	X	
	X	X			X	
			X	X	X	
		X				

WIRELESS PRO is today's most advanced electrotherapy technology offering enhanced ease-of-use and convenience while saving space and time to optimize patient treatment. It offers a new way of stimulation without constraints allowing you and your patient a total freedom of movement facilitating patient manipulation and combined exercises.

MULTI-SESSIONS Concept

Discover the MULTI-SESSIONS technology offered by WIRELESS PRO four channel devices. Four modules, two remotes and a smart docking station allow you to synchronize modules and remotes to independently manage several sessions simultaneously. Additional remotes/modules are optional.